aquarick[®]

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SAFETY PRECAUTIONS

When using electrical appliances, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER-To reduce the risk of electrocution:

- Always unplug this product immediately after using.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.

WARNING-To reduce the risk of burns, electrocution, fire, or injury to persons:

- This product should never be left unattended when plugged in.
- Close supervision is necessary when this product is used by, on, or near children or infants.
- Use this product only for its intended use as described in this manual.
Do not use attachments not recommended by the manufacturer.
- Never operate the product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center immediately for examination and repair.
- Keep the cord away from heated surfaces.
- Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, etc...
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening or hose.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Unplug this product before filling. Fill reservoir with water only unless otherwise specified by manufacturer. Do not overfill.
- The motor in aquapick has a temperature-over-heated control. When it is overheated, the power will turn 'off' automatically. After cooling down, it will resume operating normally.
- At the risk of electric shock do not remove cover(or back). Refer servicing to qualified service personnel.

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

DIRECTIONS 1

1. Unwind power cord and plug into appropriate power outlet
2. Attach jet tip by inserting it into the handle while at the same time pushing down on jet tip ring piece. Once inserted, release jet tip ring piece to lock jet tip into place.
3. Remove the water tank and open the cover to fill the tank with desired solution such as water (cool or warm) or oral rinse.
4. Align tank with mount and press gently until secure.

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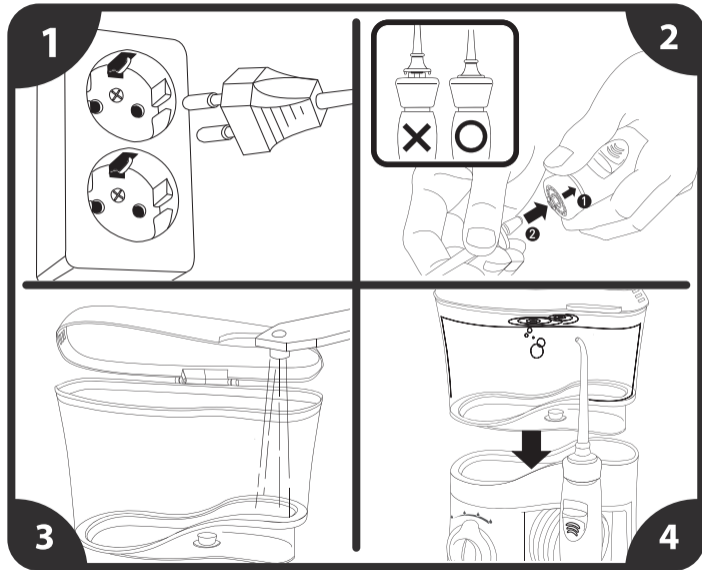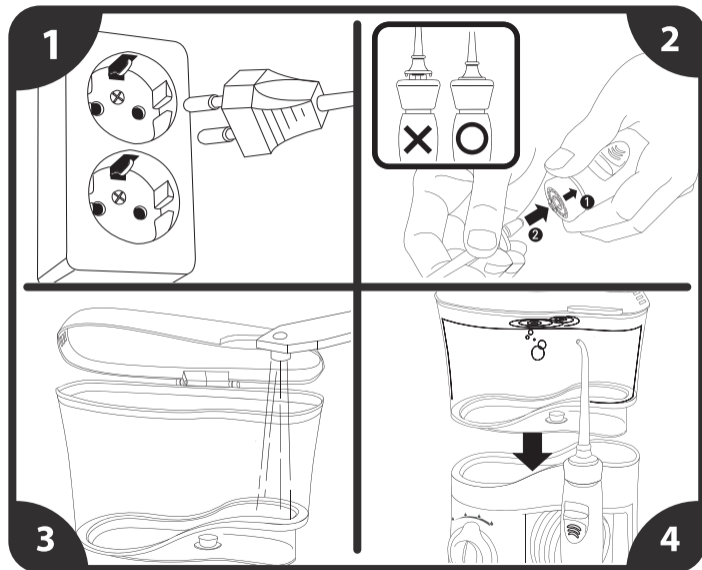DIRECTIONS 2


5. Adjust water pressure knob.
Right side water drop (strongest) / Left side water drop (weakest)
WEAK(Low) First time users/children
STRONG(High) Adjust to desired strength and comfort level.

6. Make sure the water switch button is in the off(downward) position. Press the power button on to start motor.

7. Place jet tip against desired area of mouth.
(Lean forward to allow excess water to leave mouth easily)
▪ ***Place jet tip to mouth before pushing the water switch button on.**

⚠ Notice

If there is any bleeding phenomenon after using aquapick it means that you may have periodontal disease. Continue to use for two weeks, and the bleeding should stop. If the bleeding continues, please stop usage and consult your dentist

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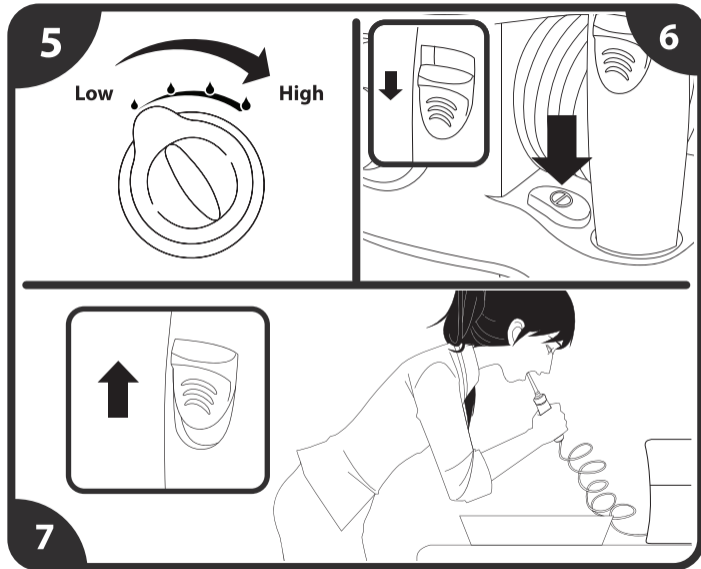
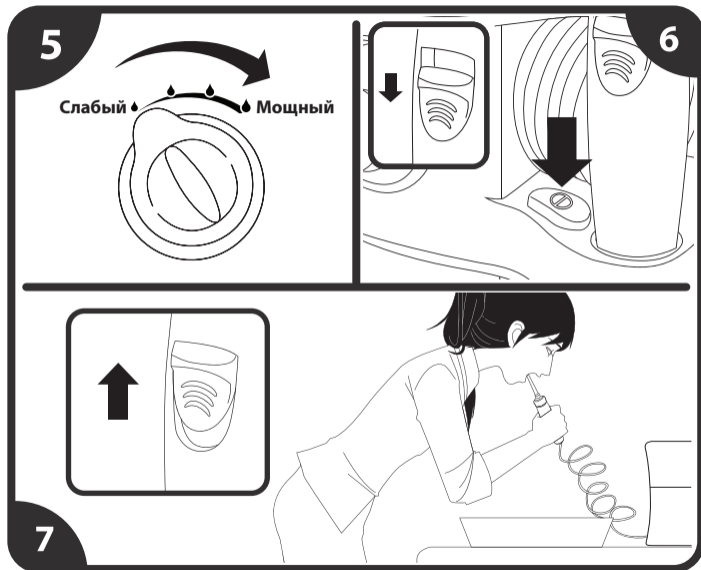
OPERATION

Washing Steps


1. Place the water tank filled with solution on the aquapick base and turn on the power button switch.
2. Slowly bend your face close over the basin, and keep your face at a 45 degree angle to the basin.
3. Press the jet tip to your mouth to begin washing.
4. Push the water switch button on handle to turn on.

How To Irrigate The Teeth & Gums

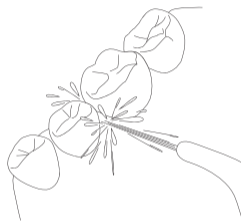1. Contact the tooth at a 90 degree angle, and move along the teeth clockwise.

2. Please cleanse sensitive area between teeth carefully.

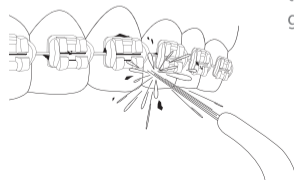3. Direct the jet stream along the gum line with the jet tip. It will flush away remaining food particles on and between teeth effectively while massaging the gums gently.

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
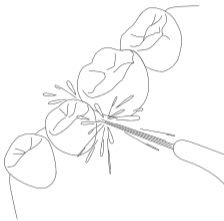
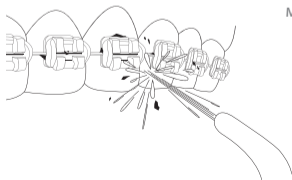
DISPOSAL 1

1. Take off the cover and pour the remaining water out.

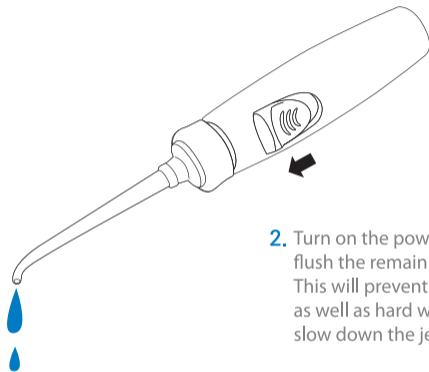2. Turn on the power for a few seconds to flush the remaining water out of the tank. This will prevent bacteria from growing as well as hard water buildup which will slow down the jet stream flow.

■ Please maintain once a week according to 'cleaning and maintenance'

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

DISPOSAL 2

- Put the handle back in the holder, Turn the power button off and pull water switch button to downward position.

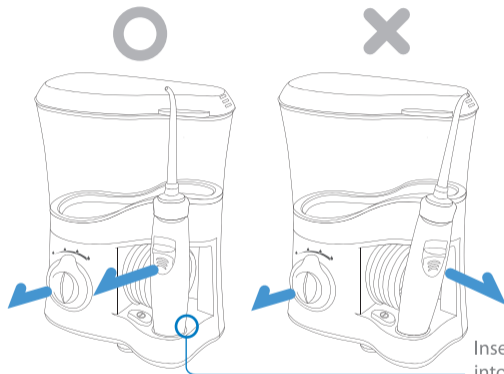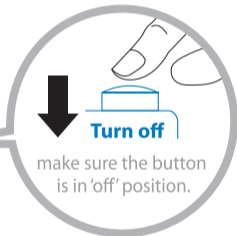- To ensure proper storage, please make sure hand piece is facing forward when inserted.



Insert the bottom of the handle into the aquapick base to ensure convenient and effective storage.

- After use, please remove jet tip and place into jet tip holders.

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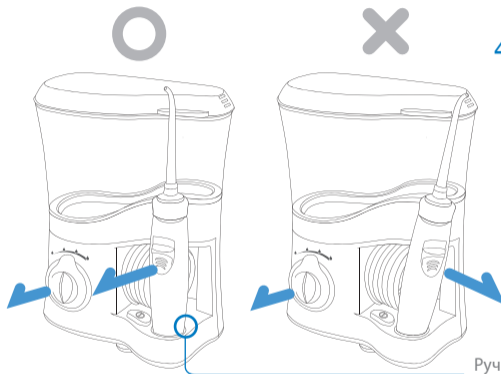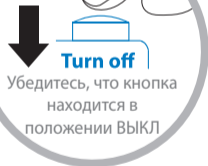TROUBLE-SHOOTING GUIDE

Please try trouble-shooting options before sending in for repair.

Problem	Cause	Corrective Action
Aquapick is not turning on	Improper connection with the power source	Unplug and plug in again
Motor is not working	Motor may be overheated causing the auto-safety system to be turned on	Wait for 20 minutes until the motor is cooled down
Pump is functioning But water is not spraying out	The water tank is empty	Refill the water tank with water
Water pressure is too low	Unsuitable water pressure	Adjust the knob to the most suitable setting <hr/> Check the water tank to see if it is placed in its correct position

Warning

Please do not take the appliance apart yourself to repair and assemble. If you want to replace a part such as a spray tube, water tank or jet tip, please contact your dealer or authorized maintenance center.

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