









۲

()

Perfect for me, ideal for my practice. With Hyperion X5 I get immediate diagnoses, I can motivate my patient and my work, and I can get the most out of my time. Just right for you.



reddot award 2015 winner

3

Hyperion X5, it's easy.

Compact and easy to use. Capable of being installed onto any wall suitable for an intraoral X-Ray. Out of our experience, this product represents the best solution for every dentist.

Hyperion X5 is the fruit of an innovative design capable of adapting to any practice, by ensuring real-time diagnosis. State-of-the-art technology, simplicity of use. This is the formula MyRay offers so as to combine image quality and speed, thereby lending value to your time and work.

Easy imaging.

- Setting-free
- Easy to use
- Real-time diagnostics







The best design for a better practice.

EASY USE

Simple and intuitive use, setting free, problem free.

FIT FOR ANY SPACE

For every intelligent space, for every requirement. One wall, zero footprints.

Easy, it fits.

A refined, clean, essential, ultra-compact design. Hyperion X5 is the smallest in the world and can be installed in any practice: all you need is one wall!

Hyperion X5 needs a minimum operational space, the same as an intraoral X-Ray. The captivating form and the elegant lines make it an agreeable element capable of adding value to your practice. Hyperion X5 does not clutter and blends into your environment. Thanks to the smart packaging solution, it can be quickly and safely installed so as to minimize installation times and costs. **Always by your side.**

•	Compact desig	n
---	---------------	---

• Wall-mounted

• Easy installation



Light and compact like an intraoral X-Ray, for a full world of possibilities. All you need is one wall.



Easy, it helps.

User-friendly. Quick. Setting-free. Hyperion X5 gathers in itself all the technological innovations that simplify the workflow of your practice. Everything proves therefore to be simple and intuitive.







Immediate access, facilitated visual control, intuitive interface. Interacting with a panoramic imager has never been so simple.

- Focus-free
- MRT (Morphology Recognition Technology)
- 3 laser guides
- Firm support with self-blocking wings
- Virtual console
- Facilitated access (even on wheelchair)


Thanks to the automatic focusing, with parameters automatically set for the patient, Hyperion X5 avoids the production of underexposed or overexposed images, thereby preserving the patient from useless radiation. With the MRT (Morphology Recognition Technology) technological innovation, no operation to set the emission parameters is in fact required. Through the convenient App*, you can select the desired type of examination in a simple and intuitive manner. Free your mind.









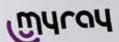

Easy, it guides.

The diagnostic features you need in just one touch. Accurate diagnoses for different requirements, from adults to children, essential programs for your needs, in a few, simple steps.

A variety of settings which suit your needs, with 18 available programs, such as: dentition with improved bitewing orthogonality, temporomandibular joints and maxillary sinuses. With the optimized orthogonal interproximal projection, you can get dentition images for a better periodontal check-up. Useful for patients unable to tolerate intraoral bitewings, as it is more comfortable and less invasive. QuickPan allows a reduction in exposure times and radiation dose, thereby safeguarding the patient's health. The MultiPAN Function offers you the possibility of choosing the panoramic imaging best suited to highlight the clinical detail of relevance from a range of different views obtained through a single examination. Ultimate flexibility.

- High quality panoramic imaging in just one touch
- 18 diagnostic programs
- Orthogonal bitewing dentition
- QuickPan (low dose)
- MultiPAN Function





Easy, it diagnoses.

Quality at the service of the clinic. The innovative technology and the automatisms inherent in the device combine so as to give birth to simply perfect panoramic images.

The high signal-noise ratio enables you to see the clinical details with utmost clarity. The broad focal layer ensures quality images in every circumstance. The filters available in iRYS make it possible to maximize the contrast and to highlight the profiles without emphasizing the noise. Thanks to an innovative filter, it is possible to optimize the perception of clinical details of diagnostic relevance by maximizing the contrast and by offering an extremely high level of detail and image clarity.

Empower your vision.


- Super Image Enhancer
- Maximum clarity
- In-depth clinical investigations









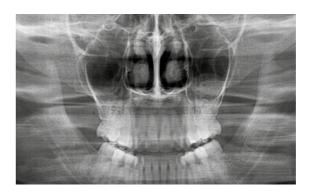

ADULT PANORAMIC IMAGING

The program of standard panoramic exposure enables a thorough or partial analysis of the patient's state by selecting the area of diagnostic relevance.

 Thorough

• Partial

CHILD PANORAMIC IMAGING

Automatically reduced exposure, same for the exposed area, thereby providing thorough diagnostic information with a restrained dose.

- Quick
- Low dose

TMJ'S (TEMPOROMANDIBULAR JOINTS)

Study of the temporomandibular state in a simpler fashion thanks to latero-lateral and posteroanterior images, 4 X-Rays with a single scan.

- Open & Closed mouth
- Sagittal & Coronal

MAXILLARY SINUSES

Characterized by a special image layer which returns X-Rays where the maxillary sinuses are clearly visible.

 Frontal

Lateral

Easy, it connects.

Managing, connecting, sharing. The image exchange and assessment has never been so quick, thanks to the Ethernet connection and the APPS for touchscreen devices.


Hyperion X5 shares its own images with all the PC's present in the practice thanks to the image management iRYS software, common to every MyRay device. It is possible to send images with a TWAIN protocol to every compatible management software. Hyperion X5 is an open system: the virtual console, available both for WINDOWS systems and as iPad APP*, enables a remote control of the machine and a visualization of the images acquired on tablet. Managing and archiving the images is practical and quick: because of that, we have conceived a powerful platform capable of interfacing with third party systems, thanks to the DICOM protocols and other communication methods. View everywhere.

REMOTE CONNECTION			

•	Ethernet connection

• iRYS Software

- Virtual console
- APPS*
- DICOM

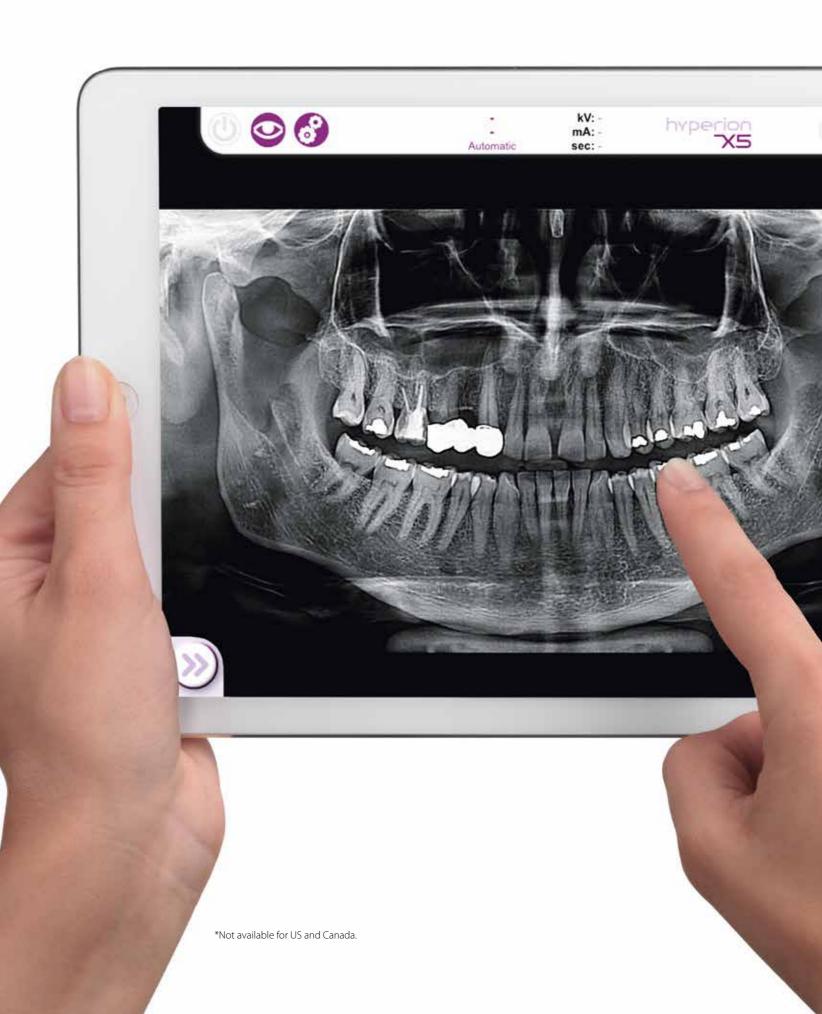

Easy, it saves.

Your time, our priority: Hyperion X5 is easy to install and is always ready for use with minimum costs and times. You, too, should choose Hyperion X5, it is simple and does not require any pre-settings.

Automatic stand-by and immediate reactivation for minimum electricity consumptions. A saving for the environment and for your budget. Thanks to the excellent quality of materials and the simplicity of the device, the assistance and maintenance needs have never been so limited, thereby ensuring maximum reliability over time. Control the state of the machine and update it even from afar with the ethernet connection. **Enjoy your time.**

- Plug&Play
- Eco friendly
- Mini budget
- Stress-free







Hyperion X5, improve your day.

Concrete benefits every day. Hyperion X5 lends value to your time, nullifies stress and maximizes your work.



It improves the

service, offering

at once an answer

to the problem by

patient's state during

interruptions. Better

workflow, a calmer

monitoring the

the treatmemt

without any

patient.

quality of the clinical



The DICOM network

ensured thanks to

protocols, available

enable the printing,

retrieval of images

and the interface

with booking lists.

with iRYS, which

archiving and

connection is



Available iPad* applications for a Wi-Fi operated remote control and quick and nimble diagnostics. Setting, start and image acquisition are within reach.

*Not available for US and Canada.






EVERYWHERE





Compact and light, Hyperion X5 calls for a minimum floor-mounted operational footprint by retaining tear strength requirements which are equivalent to an intraoral X-Ray.

Installation and calibration costs and times are restrained. The high quality of the materials calls for minimum maintenance whilst ensuring maximum reliability.



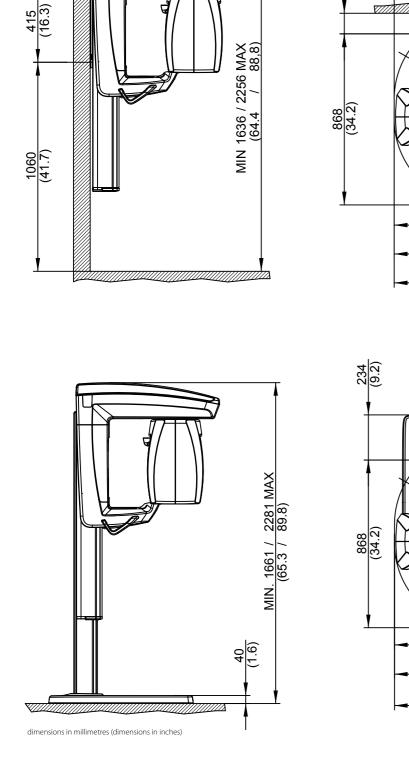

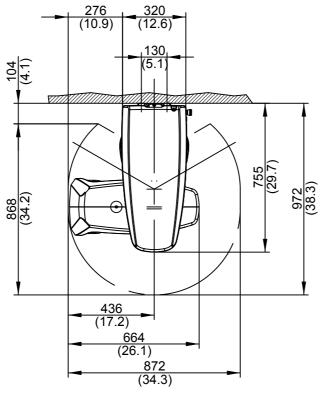

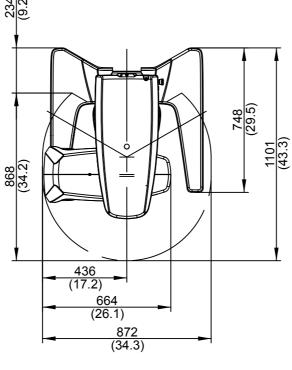

	2D IMAGES		
Туре	Adult and child panoramic imager, QuickPAN, MultiPAN, "Bitewing" Dentition, PA and LL maxillary sinuses (Right and left), Temporomandibular Joint (2 x LL +2 x PA) Open and Closed Mouth		
Patient Alignment	3 Laser Guides		
Sensor Technology	CMOS (CSI)		
Dimension	PAN standard 15 cm x 28 cm		
Image file size	max 7.5 MB		
Image resolution	5 to 7 lpmm		
Dynamic range	14 bit (16383 Gray levels)		
Enlargement	1.2 – 1.3		
Exposure time	3-12 s		
	X-RAY GENERATOR		
Type of generator	Constant potential (DC)	Constant potential (DC)	
Anode voltage	60 kV – 85 kV (step 1kV)		
Anodic current	4 mA - 15 mA		
Focal spot	0.5 mm (IEC 60336)		
Inherent filtration	> 2.5 mm Al eq. (a 85 kV)		
Exposure control	Automatic. MRT (Morphology Recognition Technology).		
Duty cycle	1:20 at full power operation (85kV, 12mA)		
	SUPPLY		
Voltage/Frequency	115 - 240 Vac, ± 10% single phase 50 / 60 Hz ± 2 Hz		
Maximum current temporary peak absorption	orary peak 20A at 115V; 12A at 240V		
Current consumption in standby mode	maximum 0.5 A		
Note	Automatic voltage and frequency adaptation		
DIMENSIONS	WALL-MOUNTED VERSION	FREE-STAND FLOOR-MOUNTED VERSION	
Required minimum	872 mm (L) x 972mm (D)	872 mm (L) x 1101 mm (D)	
Packaging dimensions (L) x (P) x (H) in mm	(1) 930x690x960 + (2) 1460x350x350	(1) 930x690x960 + (2) 1860x355x350	
Weight	70 kg (154lb)	84 kg (185lb)	
Adjustable height	2 speeds motorized		
Note	Accessible by patients on wheelchair	Self-supporting base available	
	CONNECTIVITY		
Connections	LAN / Ethernet		
Software	RYS		
Supported protocols	ported protocols DICOM 3.0, TWAIN, VDDS		
DICOM nodes	IHE certification (Print; Storage Commitment; WorkList MPPS; Query Retrieve)		

MyRay, just right for you.

Leading European company for a full range of dental imaging solutions. For all dental practices.



RXDC HyperSphere⁺ The high-technology X-Ray for maximum operational freedom.

RXDC eXTend The highly precise and versatile X-Ray with a wireless control device.



Hyperion X9 The extraoral 3-in-1 (3D, PAN, CEPH) imaging system suited to every clinic.

Hyperion X5 The compact and simple solution in the digital field of dental panoramic imaging. iRYS Best software for your 2D and 3D clinical needs.

X-pod Wireless diagnostics within everyone's reach, fitted with a touch-sensitive display.

Zen-X The fastest way of obtaining perfect intraoral images.

RXAC The solid and reliable X-Ray with AC technology for digital and film.

C-U2 High Definition (HD) videocamera with the fastest communication in the world.



www.my-ray.com

 Plant - Via Bicocca, 14/c - 40026 Imola - Bo (Italy) tel. +39 0542 653441 - fax +39 0542 653555

 Headquarters - Cefla s.c. Via Selice Provinciale, 23/a - 40026 Imola - Bo (Italy) tel. +39 0542 653111 - fax +39 0542 653344

 Cefla North America, Inc. - 6125 Harris Technology Blvd. Charlotte, NC 28269 - U.S.A. Toll Free: (+1) 800.416.3078 Fax: (+1) 704.631.4609