


# ACCEL-21

PORCELAIN FURNACE

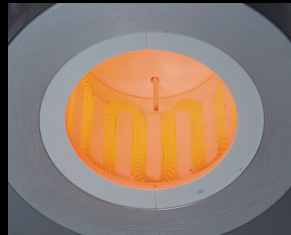

## Muffle Distributes Heat Evenly

The temperature control circuit maintains optimal firing conditions without constant user calibration.



## Rotating Firing Platform

Individual firing platforms can be manually rotated to monitor all sides of glazed or soldered pieces.



## Convenient Tool Storage

The space under the ceramic tray stand holds instruments and firing trays for easy storage.

## Various Convenient Functions

- \* Customize settings for pre-dry, drying, cooling and vacuum start temperature.
- \* Controls moisture in muffle
- \* 100 programmable memory slots

## Accel 21 Specifications

Power Requirement: AC110-120V 50/60Hz 1.4kW  
Overall Dimension: 9.8 (W) x 19.2 (H) x 14.7 (D) inch  
(250 (W) x 488 (H) x 375 (D) mm)  
Weight: 39 lbs (17.8 kg)

Vacuum Pump Outlet: AC120V, 300VA max.  
Muffle Dimension: dia. 3.6 x 2.4 inch (dia. 93 x 62mm)  
Firing Platform: dia. 3.2 inch (dia. 83mm)  
Max. Temperature: 2192°F (1200°C)