

#whdentalwerk



video.wh.com


Now at  
your dealer  
or [wh.com](http://wh.com)



# Piezo tip system.


Durable. Safe. Effective.

By professionals for  
professionals.

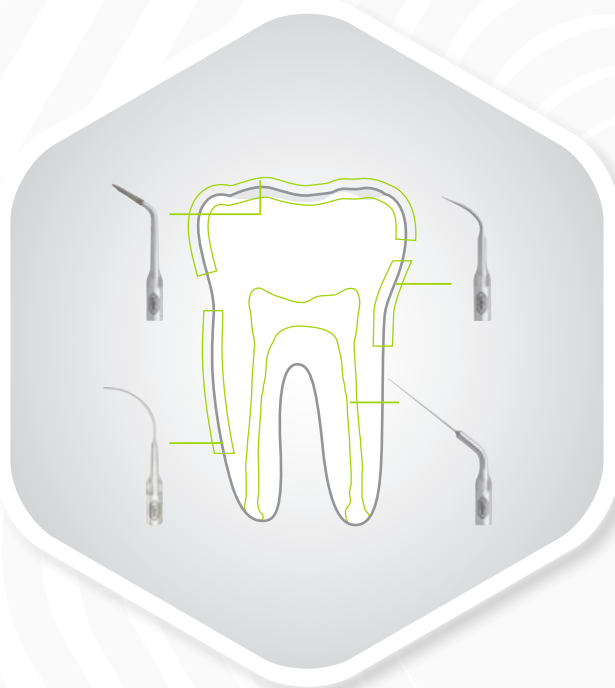

# W&H Piezo tip system

The W&H tip system has been developed in consultation with end users and specially designed for the technical requirements of W&H piezo technology. It is perfect for the requirements of the respective application areas and the individual treatment steps.

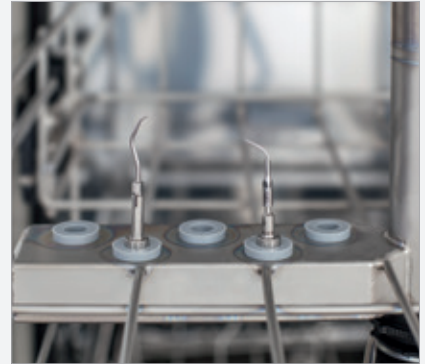

## Hygiene and Maintenance

All tips and the tip changer are thermo washer disinfectable and sterilizable.



## Fast, simple and safe handling


The W&H tip changer allows fast and safe handling. The integrated torque control prevents overtightening and leads to a safe and simple tip change.



# Universal tips

Simple, gentle and efficient supragingival calculus removal in all quadrants.

# Scaling

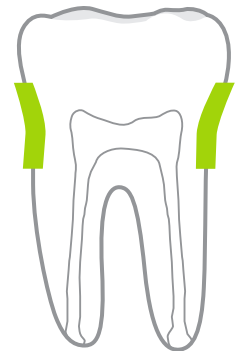

ideal 30  
max. 40

## 1UQ / 1U / 1US

For coarse depuration of supragingival calculus.

ideal 20  
max. 30

## 2UQ / 2U / 2US

For fine depuration of supragingival calculus, as well as for the treatment of recall patients.

ideal 20  
max. 30

## 3UQ / 3U / 3US

For polishing to remove supragingival calculus, particularly suitable for interdental cleaning.

ideal 30  
max. 40


## 4U

For coarse depuration of supragingival calculus, in particular for the removal of large areas of calculus.

# Periodontology tips

Paro

Its slender design facilitates optimum and also atraumatic access to the periodontal pockets.

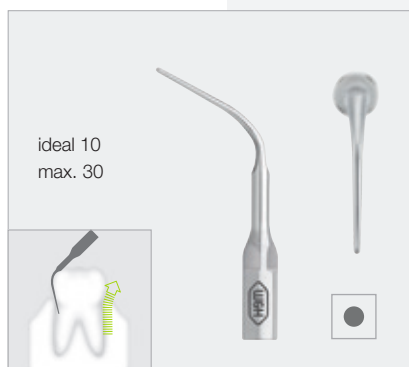


### 1PQ / 1P / 1PS

Removal of subgingival deposits, particularly suitable for deep periodontal pockets.



### 3PrQ / 3Pr

Right angled perio tip. Recommended for maintenance of moderate and deep pockets and furcations.

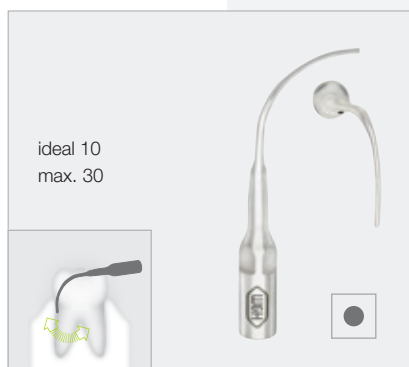

### 3PIQ / 3PI

Left angled perio tip. Recommended for maintenance of moderate and deep pockets and furcations.



### 1IQ / 1I Implant-clean



Tip for the efficient and careful cleaning of implants and the attendant restorations.

**For supportive periodontitis treatment (SPT), recall and biofilm management.**

# Periodontology tips


Paro

Diamond-coated tips for thorough cleaning of the root surface during the initial treatment of periodontitis affected patients.

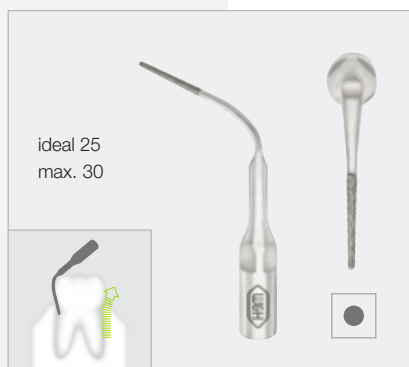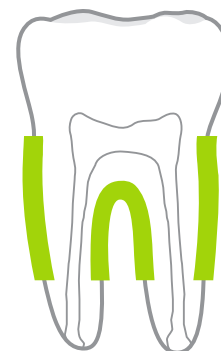


### 4PQ / 4P / 4PS

Diamonded tip for efficient cleaning of the root surface.

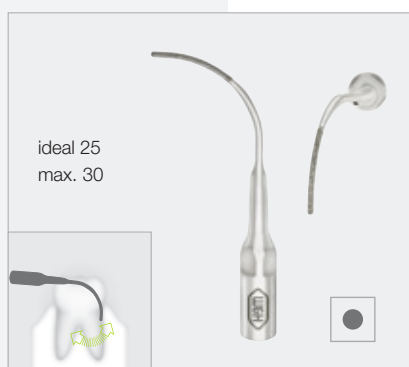


### 5PrQ / 5Pr

Right angled tip with diamond coating for efficient cleaning of the root surface.

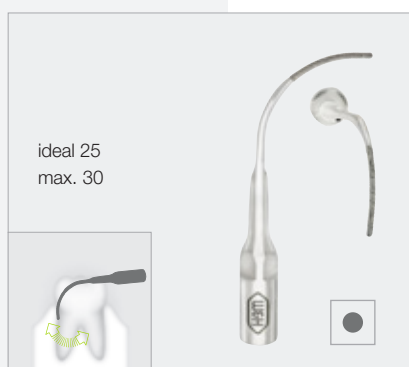

### 5PIQ / 5PI

Left angled tip with diamond coating for efficient cleaning of the root surface.

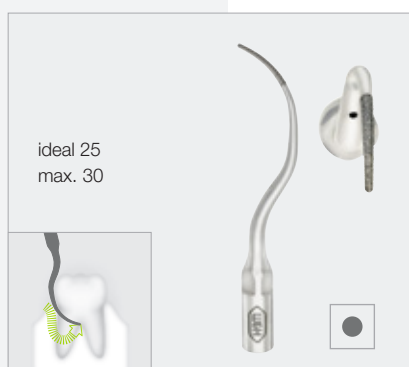


### 5PQ / 5P

Diamonded tip for efficient cleaning of furcations and concavities.


**For initial periodontitis treatment and removal of concretions.**

# Endodontic tips

Endo

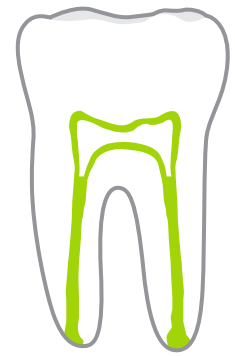

W&H offers the perfect range of endo tips whether for flushing a root canal or complex revision in an endodontic patient.

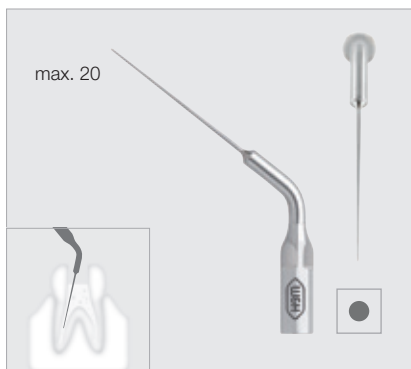

## 1E / 1ES

Non cutting instrument for the irrigation at the end of the root canal preparation.



## 2E

Tip for the lateral condensation due to thermal effect.



## 3E

Removing of broke instruments in the root canal.



## 4E

Preparation of the pulp chamber.



## 5E

Vibrating and removal of root canal posts.



## 6E

Removal of root canal fillings and root canal calcifications.


# Restoration tips

The tips are the perfect addition to your straight and contra-angle handpieces. Due to the diamond coating a high abrasion performance is achieved.

# Restoration



**1R**  
Expanding of fissures.

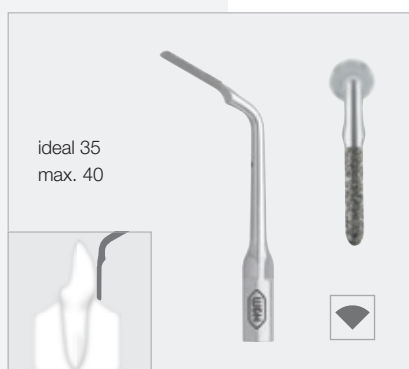

**2R**  
Diamond tip. Preparation of the veneers for clearly defined edges of preparation with maximum protection of the adjacent tooth.



**3Rd**  
Partly diamond coated tip.  
Preparation for the of the distal surface.



**3Rm**  
Partly diamond coated tip.  
Preparation for the of the mesial surface.

# Compatibility of tips with W&H handpieces

	Q-Link*	W&H/EMS/Mectron			W&H/NSK/Satelec**		
	PB-5 L Q	PB-5 L	PA-1 LED	PB-3	PB-5 L S	PB-3 LED S	PB-3 S
<b>Tigon+</b>			X			X	
<b>Proxeo Ultra</b>	X	X			X		
<b>Pyon 2</b>				X			
<b>Built-In</b>	X	X		X	X		X

\* Q-Link tips (e.g. 1UQ, 2UQ) are compatible with the PB-5 L Q.

\*\* W&H tips with Satelec thread (e.g. 1US, 2US) are compatible with the W&H handpieces PB-5 L S, PB-3 LED S and PB-3 S.

Note: EMS®, Mectron®, NSK® and Satelec® are trademarks of third parties not affiliated with W&H Dentalwerk Bürmoos GmbH. When using third-party tips, the user has to adjust the performance parameters of the W&H Piezo unit accordingly.